

## FAQ – General

*How do I select the appropriate current range setting?*

**Subject: Instrument Properties – Current Range**

### Introduction

The primary purpose of a potentiostat is to control the potential (voltage) of an electrochemical cell while measuring the current response. Current flowing through a resistor is directly proportional to the voltage difference across that resistor, as defined by Ohm's law ( $E = IR$ ). A potentiostat forces the current to flow across a resistor ( $R_m$ ) and measures current by the resulting voltage drop across  $R_m$ .


The current range defines the maximum amount of current that can be measured and the accuracy of a given potentiostat as a function of its available ranges. Princeton Applied Research offers an array of choices, all of which include the option to auto range, while the specific ranges depend on the instrument model. By selecting the current range according to best practices, the user is able to take full advantage of the potentiostat's capabilities.

Auto ranging allows the instrument to dynamically adjust the ideal range throughout the experiment. So that, at any given time, the instrument selects the lowest range for which the measured current will not exceed the range magnitude. The advantage of auto ranging is high accuracy and resolution throughout the experiment and reduced noise. One disadvantage is the possibility of data spikes, which typically arise when performing fast potential pulse/step or voltammetry experiments or when using a highly capacitive electrode. The origin of the spikes is best explained by considering how the current is measured within most Princeton Applied Research potentiostat/galvanostat systems.

An electronic circuit called a current-to-voltage converter, or simply I/E converter (Fig. 1), is responsible for measuring current. The partial block diagram in Fig. 1 shows several highly accurate resistors as part of the I/E converter, each of which represents an available current

range (in this example, 20 nA to 2 mA). Any of the available current ranges can be selected for continuous use throughout an experiment. This is referred to as a fixed current range, as opposed to auto range. In either case, the current passes across one of these resistors such that a voltage of reasonable magnitude can be measured across that resistor. Thus, currents are measured via voltage measurements across  $R_m$ . Recalling Ohm's Law, a low current (100  $\mu$ A) across a large value  $R_m$  (10 k $\Omega$ ) produces a larger voltage (1 V) compared to passing that same current across a small value  $R_m$  (10  $\Omega$ ,  $E = 1$  mV). This explains why the value of  $R_m$  increases as the current range decreases.

The "Auto" range setting uses a series of relays to switch from one  $R_m$  to another. As one current range begins to lose accuracy due to low current (thus, low voltage) across a given  $R_m$ , a relay switches to a greater magnitude  $R_m$ , thereby increasing the voltage and providing greater accuracy in measuring the lower current. For example, a

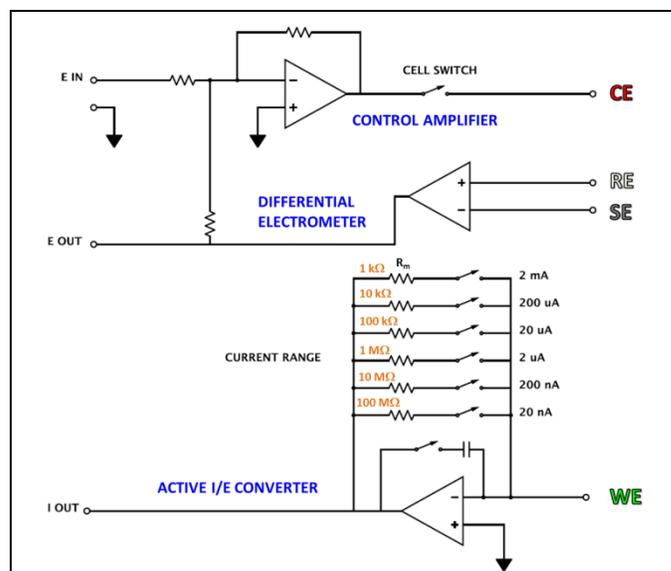

**Figure 1: Typical potentiostat block diagram**

10  $\mu\text{A}$  current flowing across the 2 mA range ( $R = 1 \text{ k}\Omega$ ) would produce a voltage of 10 mV; switching that same current across the 200  $\mu\text{A}$  range ( $R = 10 \text{ k}\Omega$ ) increases the voltage to 100 mV, which can be measured with greater accuracy and precision. When the current increases and a current overload is detected, the circuit behaves accordingly, switching to an appropriate higher range with a smaller magnitude  $R_m$ .

## Practical Examples


With a basic understanding of the potentiostat's hardware that is responsible for measuring current, we can use Cyclic Voltammetry (CV) to demonstrate how the Current Range setting can affect the resulting data. VersaStudio allows the user to select between two voltammetry Actions: Staircase Voltammetry and Cyclic Voltammetry. The data included in this document were collected using the latter Action, which is designed to mimic an analog ramp generator by automatically making the potential step size as small as possible, using a maximum of 1000 points per scan or 2000 per CV. Note that Staircase Voltammetry allows the user to control the number of points collected by defining a step height and time. All CVs were collected using a micro-cell kit (K0264), Au milli electrode (G0227) and Pt wire counter electrode (K0266). All potentials reference a Ag/AgCl reference electrode. The solution was 1 mM  $\text{K}_3[\text{Fe}(\text{CN})_6]$  in 0.1 M KCl.

Fig. 2 shows 3 CVs run on the same cell at 100 mV/s, using various settings. Fig. 2 (A) was collected using the I Filter set to 10 Hz and Current Range set to "Auto", while Fig. 2 (B) used the same I Filter and a 20  $\mu\text{A}$  fixed current range. The I Filter and Current Range were both set to "Auto" in Fig. 2 (C), where the "Auto" I Filter setting is 1 kHz. Fig. 3 shows the *Experiment Properties* window with the *Advanced Properties* view enabled, where these settings were assigned.


Let's first consider why there are spikes in the data shown in Fig. 2 (A). If the raw data are examined more closely, it is clear that the spikes occur when the current range (*I Range*) switches from one range to the next.

Fig. 4 shows the current spike at  $\sim 270 \text{ mV}$  observed in the negative going scan from Fig. 2 (A) and the corresponding raw data. Prior to the spike, the current is measured on the 2  $\mu\text{A}$  range as  $-1.833 \mu\text{A}$ . When the current increased beyond 2  $\mu\text{A}$ , the relay switched to the 20  $\mu\text{A}$  range, resulting in a spike to  $-15.865 \mu\text{A}$ . **The**

(A) I Filter = 10 Hz, Current Range = "Auto"

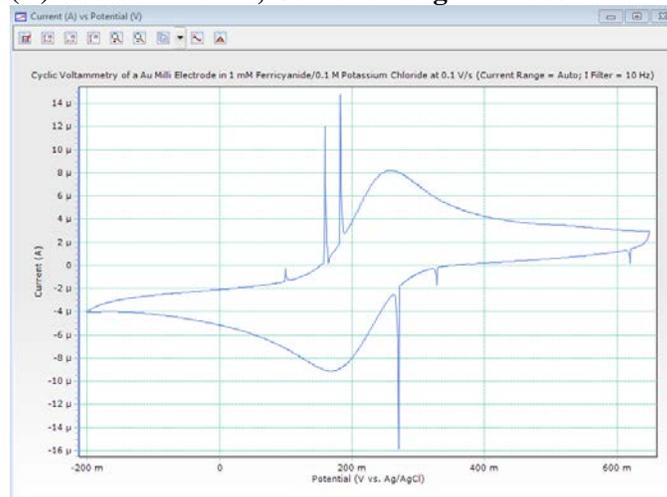

(B) I Filter = 10 Hz, Current Range = 20  $\mu\text{A}$

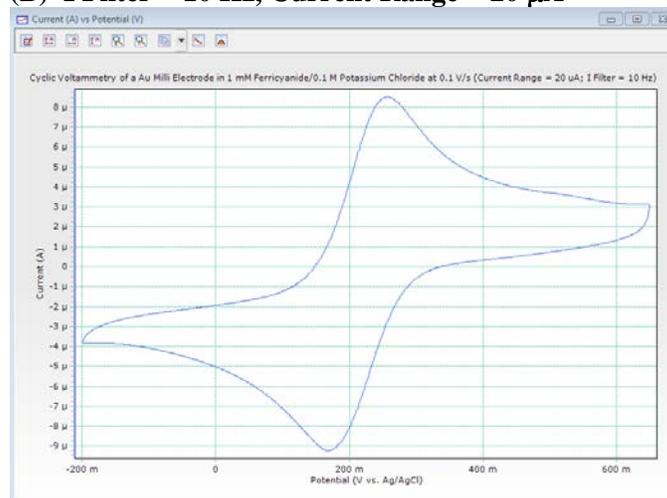

(C) I Filter = 1 kHz ("Auto"), Current Range = "Auto"


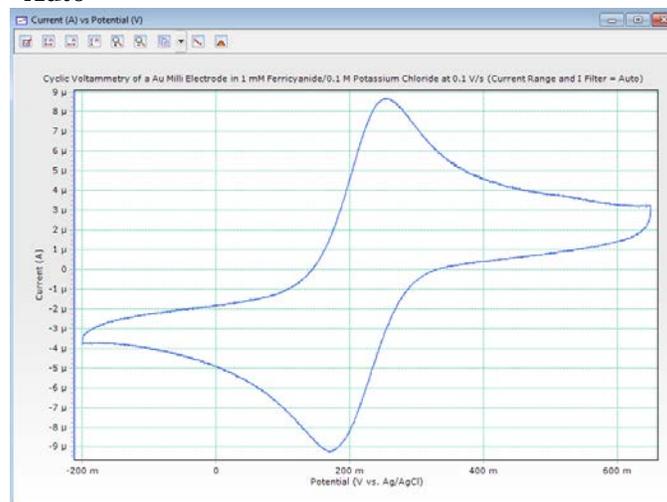


Figure 2: CVs of Au in 1 mM  $\text{K}_3[\text{Fe}(\text{CN})_6]$  in 0.1 M KCl at  $v = 0.1 \text{ V/s}$

**spike is a measurement artifact, and is not a spike that occurred at the cell.**

This is demonstrated in Fig. 5, which is a plot of E versus time for the same CV in Fig. 2 (A) and Fig. 4. Clearly there are no potential spikes corresponding with the current spikes, which verify that the cell was not impacted. Right before the range change,  $-1.833 \mu\text{A}$  was measured across a  $1 \text{ M}\Omega$  resistor, resulting in  $-1.833 \text{ mV}$ . The range was quickly changed via fast electronic relays and the next data point should have produced about  $-0.2 \text{ mV}$  for  $\sim 2 \mu\text{A}$  flowing across a  $100 \text{ k}\Omega$  resistor ( $20 \mu\text{A}$  range); however, the voltage takes time to adjust from the previous  $-1.8 \text{ mV}$  to the new  $-0.2 \text{ mV}$ . Any data collected during this transition will be obscured and should be considered measurement artifacts. The time taken to regain stability is indicative of the time required for the measured response (i.e. voltage drop across  $R_m$ ) to pass through the circuitry, which depends on the combined time constants of the circuit. By imposing a  $10 \text{ Hz}$  I Filter, we have increased the time constant relative to that imposed by the “Auto” setting ( $1 \text{ kHz}$  filter). In this case it takes about 7 data points or  $60 \text{ ms}$  to regain stability and any data points taken during this transition are invalid.

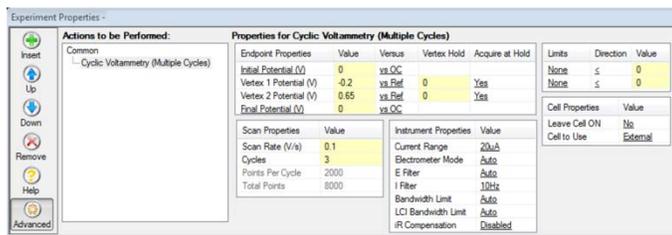


**Figure 3: Experiment Properties window**

Notice that the spike in Fig. 4 is in the direction of increasing current. This directionality results because the transition is from a smaller to a larger current range (higher value  $R_m$  to lower value  $R_m$ ). On the  $2 \mu\text{A}$  range,  $-1.833 \mu\text{A}$  produces  $-1.833 \text{ V}$ . At the instant the range changes,  $-1.833 \text{ V}$  appears across  $R_m$  ( $100 \text{ k}\Omega$ ) for the  $20 \mu\text{A}$  range producing an inaccurately large current reading ( $-15.865 \mu\text{A}$ ). Thus, the direction of the “spike” (either larger or smaller than the true reading) depends on whether the range change is a range-up or a range-down during the auto range process.

Fig. 2 (B) shows a CV of the same cell collected using a fixed  $20 \mu\text{A}$  range. There are no spikes since the current range does not change. Spikes are also absent from Fig. 2 (C) which was collected using auto current range, but with a faster  $1 \text{ kHz}$  I Filter. The faster filter results in a

shorter time constant, meaning the measured voltage passes through the circuitry more rapidly.

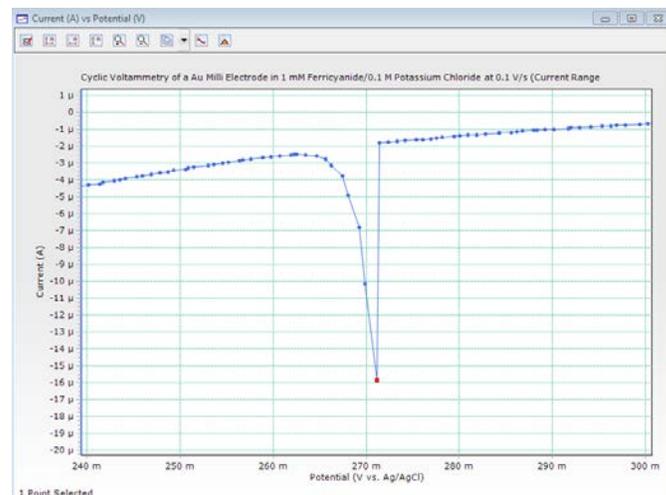


Data View - 1999 Points (Segments 2)

Potential (V)	Current (A)	Elapsed Time (s)	I Range (A)
276.955 mV	-1.558 $\mu\text{A}$	34.598752	2 $\mu\text{A}$
276.035 mV	-1.602 $\mu\text{A}$	34.607256	2 $\mu\text{A}$
275.421 mV	-1.645 $\mu\text{A}$	34.61576	2 $\mu\text{A}$
274.195 mV	-1.692 $\mu\text{A}$	34.624264	2 $\mu\text{A}$
273.275 mV	-1.737 $\mu\text{A}$	34.632768	2 $\mu\text{A}$
272.354 mV	-1.784 $\mu\text{A}$	34.641272	2 $\mu\text{A}$
271.434 mV	-1.833 $\mu\text{A}$	34.649776	2 $\mu\text{A}$
271.128 mV	-15.865 $\mu\text{A}$	34.65828	20 $\mu\text{A}$
269.901 mV	-10.159 $\mu\text{A}$	34.666784	20 $\mu\text{A}$
269.287 mV	-6.824 $\mu\text{A}$	34.675288	20 $\mu\text{A}$
268.061 mV	-4.882 $\mu\text{A}$	34.683792	20 $\mu\text{A}$
267.447 mV	-3.765 $\mu\text{A}$	34.692296	20 $\mu\text{A}$
266.22 mV	-3.13 $\mu\text{A}$	34.7008	20 $\mu\text{A}$
265.607 mV	-2.778 $\mu\text{A}$	34.709304	20 $\mu\text{A}$
264.687 mV	-2.595 $\mu\text{A}$	34.717808	20 $\mu\text{A}$
263.46 mV	-2.511 $\mu\text{A}$	34.726312	20 $\mu\text{A}$
262.54 mV	-2.486 $\mu\text{A}$	34.734816	20 $\mu\text{A}$
262.233 mV	-2.499 $\mu\text{A}$	34.74332	20 $\mu\text{A}$
261.926 mV	-2.531 $\mu\text{A}$	34.751824	20 $\mu\text{A}$
260.7 mV	-2.575 $\mu\text{A}$	34.760328	20 $\mu\text{A}$
259.779 mV	-2.627 $\mu\text{A}$	34.768832	20 $\mu\text{A}$

**Figure 4: Close-up view of a spike and the corresponding raw data from CVs plotted in Fig. 2 (A).**

The combination of time constants and the sampling rate determine if current spikes are evident in the data. Fig. 6 shows CVs of the same cell taken at  $2 \text{ V/s}$  with the  $1 \text{ kHz}$  I Filter, the same filter used in Fig. 2 (C). Current spikes and large sections of overloads (yellow highlighted data) are evident in the CV taken on “Auto” range (blue scan),

which is in contrast to Fig. 2 (C). Since the Cyclic Voltammetry Action collects 2000 points per cycle, the sampling rate (points/s) at 0.1 V/s (Fig. 2 (C)) was 20 times slower than that at 2 V/s. When auto ranging, the faster the sampling rate, the greater the chance that a spike will be captured in the data. As expected, the CV taken with a fixed 200  $\mu\text{A}$  range (red scan) is smooth. These results emphasize that there are time constants inherent to the measurement circuitry even without using the slower 10 Hz I Filter. These are due to the capacitance of the cell combined with the resistance of  $R_m$ , for which the time constant increases for larger values (lower ranges). Another factor is the relays themselves, which take milli-seconds to switch. If data are being acquired at micro-second speeds, spikes can occur when auto range is enabled.

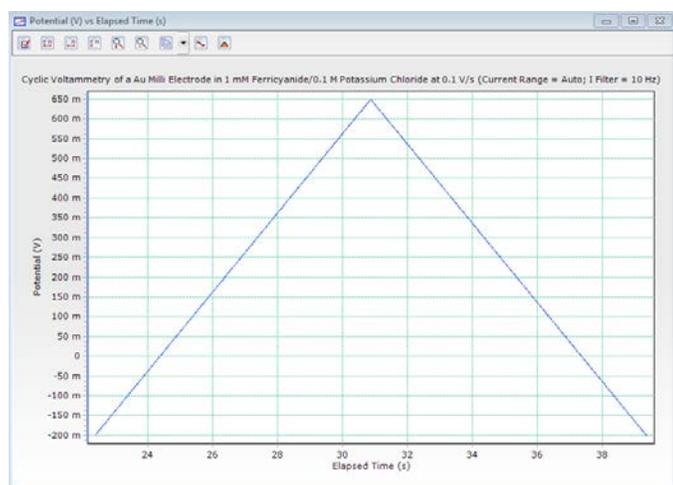


**Figure 5: Potential versus time plot for CVs collected on Auto range (Fig. 2 (A));  $v = 2$  V/s, I Filter = Auto**

## Further Considerations

The VersaStudio Manual recommends using a fixed current range if performing fast pulse/step experiments or scans greater than 20 mV/s. This general suggestion, as indicated by the data presented above, is not a hard-fast rule. The maximum scan and pulse/step rates for which auto ranging produces high-quality data is dependent upon the experimental conditions, including the cell's capacitance. Higher capacitance cells are more likely to result in data spikes. The best practice is to use "Auto" range whenever warranted. If spikes and/or overloads occur that disrupt accurate data analysis and/or when "presentation quality" results are required, repeat the experiment on a fixed range, determined by the initial auto ranging experiment, or sample data at a slower rate. For certain high capacitance samples (thick oxide on the

surface of an alloy) spikes will still occur even when using a slow rate. In this situation, it is recommended that the spikes be removed from the data set. This can be accomplished by going to *Data/Delete.../Delete Selected Only or Delete All Overloaded*.

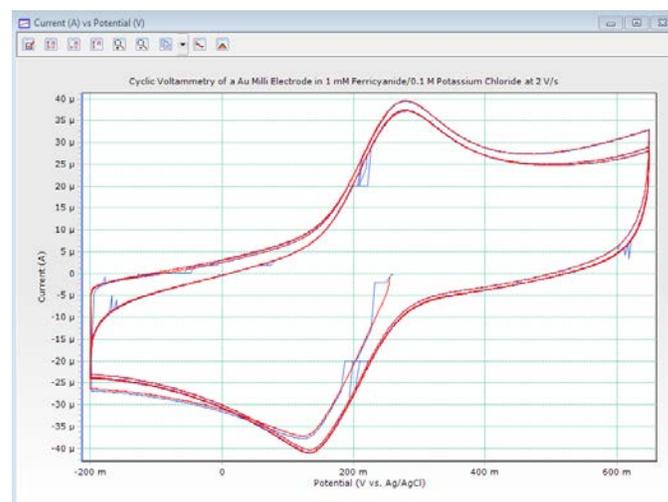

Data View - 7826 Points (All Segments)

Potential (V)	Current (A)	Elapsed Time (s)	I Range (A)	Comment
260.086 mV	-60.286 nA	0.00023	200 nA	
260.086 mV	-163.143 nA	0.00046	200 nA	Overload
259.779 mV	-201.027 nA	0.00069	200 nA	Overload
259.779 mV	-201.027 nA	0.00092	200 nA	Overload
259.166 mV	-201.027 nA	0.00115	200 nA	Overload
260.086 mV	-201.027 nA	0.00138	200 nA	Overload
258.859 mV	-201.027 nA	0.00161	200 nA	Overload
258.859 mV	-201.027 nA	0.00184	200 nA	Overload
258.553 mV	-201.027 nA	0.00207	200 nA	Overload
259.473 mV	-201.027 nA	0.0023	200 nA	Overload
259.166 mV	-201.027 nA	0.00253	200 nA	Overload
258.553 mV	-201.027 nA	0.00276	200 nA	Overload
258.859 mV	-201.027 nA	0.00299	200 nA	Overload
258.246 mV	-201.027 nA	0.00322	200 nA	Overload
258.246 mV	-201.027 nA	0.00345	200 nA	Overload
257.939 mV	-201.027 nA	0.00368	200 nA	Overload
258.553 mV	-201.027 nA	0.00391	200 nA	Overload
257.939 mV	-201.027 nA	0.00414	200 nA	Overload
258.246 mV	-201.027 nA	0.00437	200 nA	Overload
257.633 mV	-201.027 nA	0.0046	200 nA	Overload
249.045 mV	-2.01 $\mu\text{A}$	0.00575	2 $\mu\text{A}$	Overload
248.125 mV	-2.01 $\mu\text{A}$	0.00598	2 $\mu\text{A}$	Overload
247.511 mV	-2.01 $\mu\text{A}$	0.00621	2 $\mu\text{A}$	Overload
246.591 mV	-2.01 $\mu\text{A}$	0.00644	2 $\mu\text{A}$	Overload
247.205 mV	-2.01 $\mu\text{A}$	0.00667	2 $\mu\text{A}$	Overload
247.205 mV	-2.01 $\mu\text{A}$	0.0069	2 $\mu\text{A}$	Overload
246.591 mV	-2.01 $\mu\text{A}$	0.00713	2 $\mu\text{A}$	Overload
246.284 mV	-2.01 $\mu\text{A}$	0.00736	2 $\mu\text{A}$	Overload
245.364 mV	-2.01 $\mu\text{A}$	0.00759	2 $\mu\text{A}$	Overload
244.444 mV	-2.01 $\mu\text{A}$	0.00782	2 $\mu\text{A}$	Overload
244.751 mV	-2.01 $\mu\text{A}$	0.00805	2 $\mu\text{A}$	Overload
244.137 mV	-2.01 $\mu\text{A}$	0.00828	2 $\mu\text{A}$	Overload

**Figure 6: CVs collected on Auto range (blue) and a fixed 200  $\mu\text{A}$  range (red);  $v = 2$  V/s, I Filter = Auto**

There are some cases where auto ranging should always be avoided. These include very fast scans ( $>1$  V/s, i.e. Fig. 6, blue scan) and pulsing or double-step techniques (such as Square Wave Voltammetry or charge measurement Actions), since major gaps in data will occur due to time losses associated with the range change itself.

Special caution should be taken when performing experiments with low input impedance cells such as batteries, fuel cells or super capacitors. A low impedance cell is susceptible to fast changes in current, orders of magnitude in size, with only small changes in potential. The potentiostat is measuring this current and, when set to "Auto" range, attempting to use the appropriate range. Large, fast current changes can lead the hardware to fluctuate between current ranges, causing excessive overloads that produce unnecessary wear on the mechanical and electronic components of the measurement circuit. Therefore, it is recommended that the "Auto" setting be avoided if possible and that all experiments be conducted using the highest current range (2A, for example) first, only switching to lower ranges if deemed appropriate for the experiment. This is a situation where the Auto Current Range Setup Action under the *Advanced Actions* tab can be used to set an initial, minimum and maximum current range. This is achieved by inserting this Advanced Action as a step prior to your experiment Action. A detailed explanation is provided in the VersaStudio manual.