


## Advanced Instrumentation for Bioimpedance Measurements

A. J. Hinton and B. Sayers:  
Solartron, Victoria Rd., Farnborough, Hampshire UK

Bioimpedance is rapidly gaining popularity in a wide field of bio research applications including investigations into skin hydration, dental decay, body fat content, tissue ischemia, electrode / gel performance, food freshness, plant / tree growth.....

Bioimpedance is a non-destructive technique which makes use of low level electrical signals such as those which are already present in bio-materials. Accurate measurement<sup>(1)</sup> of the impedance of bio-materials over a broad frequency range yields valuable information about the electrical properties of the material. Analysis of these measurements can provide valuable diagnostic information which may be used, for example, for checking organ integrity during transplants or for tumour investigations. In addition, the study of frequency related phenomena as a function of temperature<sup>(2)</sup> is regularly used in the laboratory testing of bio materials.

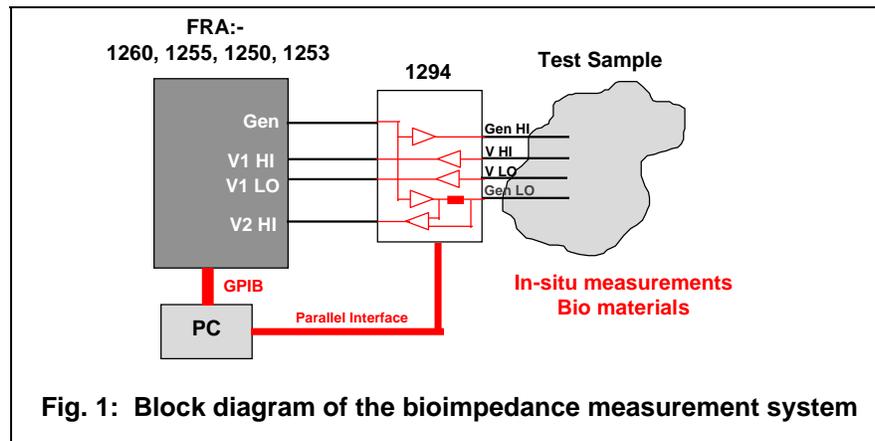

Research studies are frequently performed directly on live subjects (*in vivo* testing), for example the study of skin impedance to determine tissue healing rates or to examine the effects of creams and lotions. It is important for *in vivo* studies that impedance test equipment is designed to comply with safety standards such as IEC601. This standard outlines safe current levels which must not be exceeded, safety of equipment in the event of component failures and earth leakage requirements.

Some examples of bio-impedance applications are listed below:-

- dental research, used for the detection of decayed or cracked enamel<sup>(3, 4)</sup>.
- assessing the extent of ischemia in organ transplants; ischemia is a process by which the organ undergoes progressive change after isolation from the blood supply. Eventually, a point is reached where resuscitation is no longer possible. Impedance techniques are used to monitor the extent of tissue damage, organ integrity and the possibility of successfully reviving the organ<sup>(5, 6)</sup>.
- detection and study of tumours; the treatment of tumours via loco-regional hypothermia leads to changes in the tissue impedance due to shifting in the charged fluids present. Minimal invasive procedures such as impedance are vital in tumour assessment<sup>(7-9)</sup>.

### Other Bio-impedance Applications

- in-vivo muscle and tissue studies<sup>(10-11)</sup>
- electrode skin studies (tissue healing rates, etc.)<sup>(12-15)</sup>
- dermatological applications (medical creams / lotions)<sup>(16)</sup>
- drug delivery rates
- pacemaker development
- blood cell analysis, monitoring of viral effects on cell structure<sup>(17)</sup>
- biotechnology research, food, pharmaceuticals<sup>(18)</sup>

## Measurement Techniques

Bioimpedance measurements, whether *in vivo* or *in vitro*, impose difficult requirements on the instrumentation, including<sup>(19,20)</sup> the necessity for wide operational frequency range, accurate four terminal measurement performance, safety when connecting to live subjects and the elimination of stray impedances due to cables and electrodes. In addition, some measurements are made even more difficult due to the non-uniformity and non-linearity of the materials themselves.

The specification of the Solartron 1294 Impedance Interface is achieved by the use of various techniques including 4-terminal connections to the sample, balanced generator and driven shield voltage measurement connections. In addition, the 1294 makes use of a sensitive multi-range current to voltage converter which allows the measurement of very low current levels which are commonly experienced in high impedance analysis or when using micro-electrodes.

The 1294 operates in conjunction with a frequency response analyzer (FRA), which provides the AC stimulus signal and correlation analysis of the output signals from the 1294. The use of correlation analysis is very important for the rejection of harmonics which result from non-linearities in the samples being measured.

The combination of instruments and associated software provides:-

- The ability to test in-vivo or in-vitro samples
- Impedance measurement range to >100 Gohm (2-terminal mode), 100Mohm (4-terminal mode)
- Wide frequency range 10 $\mu$ Hz to >1MHz (IEC601 current limited output 100Hz - >1MHz)
- Increased measurement accuracy by the use of 4-terminal driven shield connections and balanced generator techniques
- Temperature control using external controller (e.g. LakeShore 340 / Oxford Instruments) and cryostat (e.g. Oxford Instruments)
- Flexible PC software allowing the control of complex experiments at a range of frequencies, stimulus levels and temperatures, and the ability to plot results in a wide variety of formats.

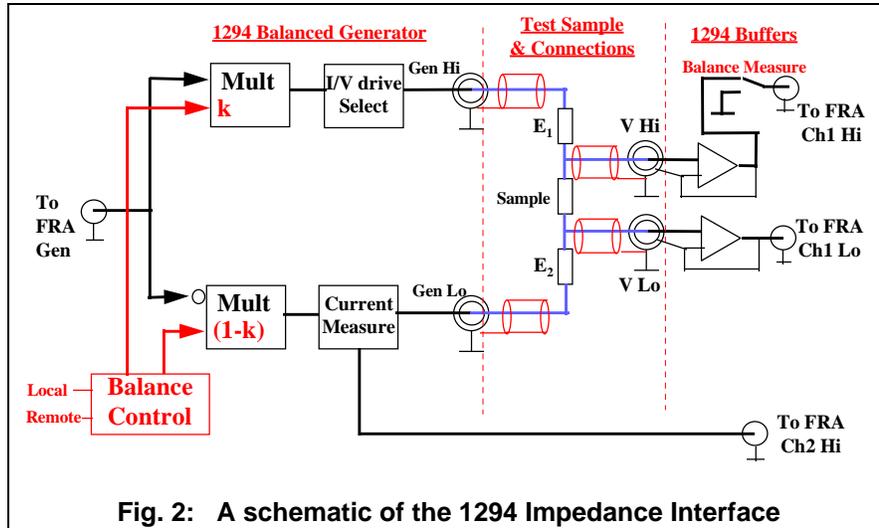


Fig. 2: A schematic of the 1294 Impedance Interface


The schematic (Fig. 2) shows some of the features which have been incorporated to provide accurate and reliable measurement of electrical impedance. The sample is positioned between V Hi and V Lo, either side are electrode impedances represented by E<sub>1</sub>, E<sub>2</sub>. These impedances are due to the sample connections which are present in many bioimpedance applications. There are also other stray impedances due to connections and the use of hydro gels (medical research), which may introduce imbalance in the electrode impedances which need to be considered and compensated.

The main features included in the 1294 are:-

- low capacitance voltage measurement inputs using driven shields
- balanced generator for better common mode rejection
- built in attenuator for low voltage / current outputs
- accurate and sensitive current to voltage converter
- bio protected sample connections conforming to IEC601

### Measurement techniques (4-terminal connections)

For bioimpedance measurements, the electrodes which connect the instrumentation to the sample may be high impedance compared to the impedance of the sample itself. If conventional 2-terminal measurement techniques are used, there is no way to measure the impedance of the sample without including the impedance of the electrodes. This can lead to inaccuracies in the measurements, additional “features” on the impedance plots and to problems interpreting the data.

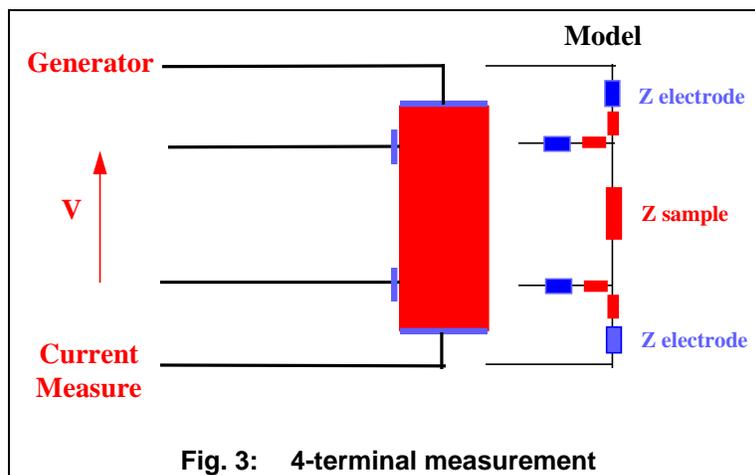

The four electrode connection technique (fig. 3) uses separate electrodes for current stimulus and voltage measurement. This allows measurements of the sample on its own and minimises errors which would otherwise be introduced by the impedance of the current carrying electrodes. If 4-terminal techniques are employed, no current flows through the voltage measurement electrodes, which means that no voltage is dropped across these electrodes, leading to more accurate voltage measurements across the sample.

In addition, there may be localised disturbances where the current is injected into the sample (GenHi and GenLo). Four terminal techniques allow the voltage measurement electrodes to be placed well away from these disturbances in a region of linear electric field giving increased measurement accuracy.

### Driven shield connections to the sample

Accurate measurements at high frequency are a specific problem due to errors introduced by input and cable capacitance. One technique which has been widely used to reduce the effects of capacitance is to position electrometer buffer amplifiers close to the sample. These buffers are able to drive the cable and input capacitance of the equipment and therefore reduce errors in voltage measurements. However, these external amplifiers require a power supply which usually involves extra cabling. In addition, where the temperature of the sample is to be varied, the accuracy of measurements may be effected since the buffers must be positioned close to the sample and are therefore subject to the same temperature variations.

An alternative solution is the use of driven shield cables. This technique replicates the signal waveform (which appears on the cable inner), onto the cable shield in order to minimise leakage current flow between the cable inner and the shield. Since no current flows between the cable inner and shield, the impedance appears to be very large and therefore the effects of the cable and input capacitance are minimised.

This method allows the high impedance buffers to be kept within the instrument. As an external power supply unit is not required for the electrometers, the cabling is reduced and the sample temperature may be changed without affecting the accuracy of the results.

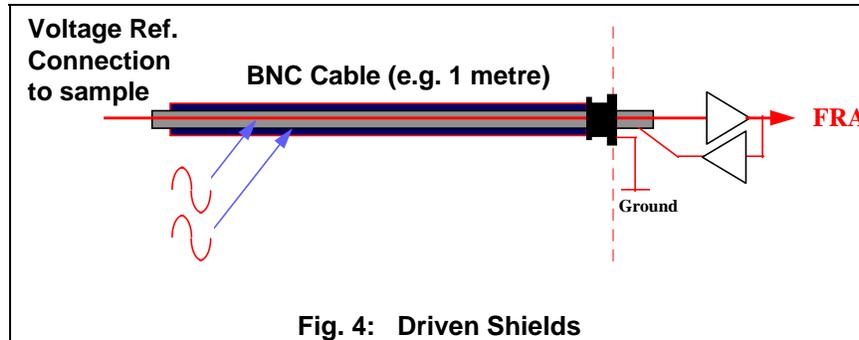


Fig. 4: Driven Shields

### Balanced Generator

One problem associated with four terminal measurements on bio-materials is the requirement to study the impedance of the sample in the presence of relatively high electrode impedances (i.e. to reject the voltages across the electrodes in order to obtain accurate measurements of the voltage across the sample itself).

Figure 5 shows a typical measurement situation where the electrode impedance is ten times higher than the sample impedance which is required to be measured. In the example shown, the sample impedance is simulated by a 1Kohm resistor, and the electrode impedances are simulated by two 10kohm resistors. Using conventional measurement techniques, the AC stimulus voltage is applied to GenHi while GenLo is grounded, and the current through the sample and voltage drop across the sample (between V Hi and V Lo) are measured in order to compute the impedance of the sample.

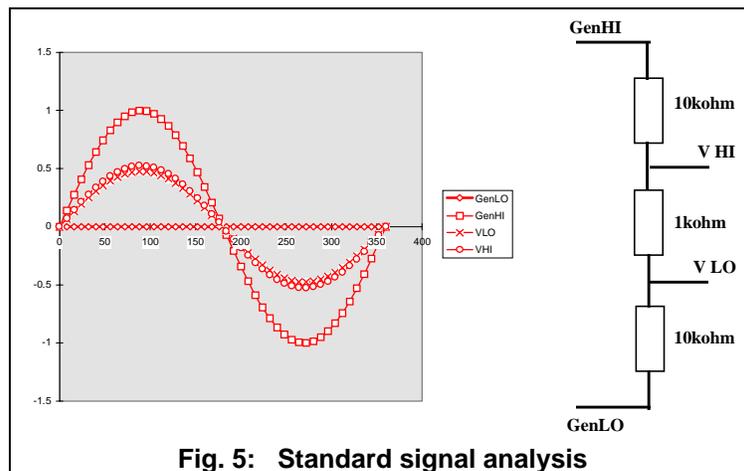
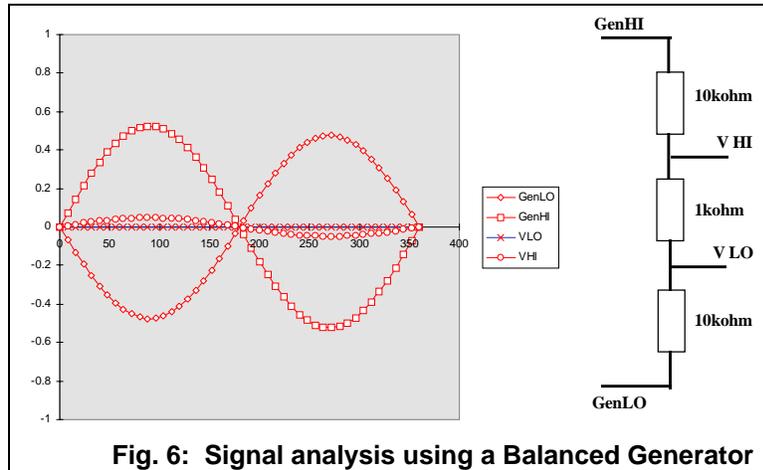


Fig. 5: Standard signal analysis

In this example, the FRA is required to measure a small voltage difference between two relatively high voltage signals which appear at V Hi and V Lo. For example, if 1 volt is applied between GenHi and GenLo, less than 50 mVolts appears across the sample (V Hi - V Lo) whereas the voltage on the V Lo measurement connection is approximately 500mV. The FRA is therefore required to measure a 50 mVolt “difference” signal in the presence of 500 mVolts which leads to errors which are referred to as “common mode” errors. In addition, the relatively high voltage on the V Lo connection causes some current which has passed through the sample impedance to leak to earth via the input / cable capacitance on V Lo instead of being measured by the current measurement circuit.

The use of a balanced generator (see Fig. 6) reduces errors due to common mode voltages and earth leakage allowing more precise measurements of the sample impedance. This is achieved by adjusting the GenHi and GenLo signals in order to provide a balanced stimulus to the sample which has the effect of making the voltage which appears at V Lo as close as possible to earth voltage, (i.e. zero volts).



**Fig. 6: Signal analysis using a Balanced Generator**

The balanced generator can also cope with extremely difficult measurement situations where for instance the electrode impedances are not equal. This is also typical for bio-impedance measurements where it is difficult to obtain reproducible electrode contacts onto skin. In this case the signals on GenHi and GenLo are set to different voltage levels in order to again achieve zero volts on the V Lo connection.

**Software capabilities**

The software developed for use with this system is compatible with Windows 3.1 / 95/ NT4. It contains a variety of functions which allow maximum flexibility in terms of measurement and presentation of results, including:-

- Control of complex experiments involving multiple frequencies, stimulus levels and temperatures
- Presentation of results in a wide variety of formats including impedance, admittance, capacitance; with the ability to plot these as a function of frequency, stimulus level or temperature
- Multiple graph overlays allowing easy comparison of results taken at different frequencies or temperatures; or overlays of previously collected data
- Driver for Oxford Instruments and LakeShore Temperature controllers

**The future of bio-impedance measurements**

Electrical Impedance Spectroscopy is a non-invasive, non-destructive technique which gives a great deal of information about the properties of the materials under investigation. The new generation of bio-impedance analysis equipment with its advances in measurement technology allows more accurate four terminal measurements of a wide range of bio-materials to be obtained. The adverse effects of stray impedances and electrode impedance are minimised by the use of driven shield and balanced generator techniques.

The technology incorporated into the 1294 Impedance Interface will push the boundaries of measurement capabilities for biomaterial characterisation. The increasing range of applications provides an indication of the vital data that this instrumentation can provide in bio-materials analysis. The trend in biomedical and clinical research is to continue to open up further areas which will benefit from electrical characterisation techniques.

## References

1. Hinton A. J and Sayers B, Impedance measurement: Sine Correlation, Solartron Tech File No. 4, 1998.
2. Gersing E, Kruger W, Osypka M and Vaupel P, Problems involved in temperature measurements using EIT, *Physiol. Meas.*, **16**, A153 - 160, 1995
3. Vandernoot T J and Levinkind M, Ac impedance characteristics of human dental enamel and dentine, *J. Electroanal. Chem.*, **300**, 191-198, 1991
4. Levinkind M, Vandernoot, T J and Elliot J C, Electrochemical Impedance Characterisation of Human and Bovine Enamel., *J. Dental Res.*, 1806 - 1811, 1990
5. Gersing E, Bach F, Gebhard M M, Kehrer G, Meissner A, Bretschneider H J, Monitoring alterations caused by ischemia in organ tissue by electrical impedance spectroscopy., NSC-BME, Antwerp, Belgium, 1990.
6. Gersing E, Burger E, Gebhard M M, Kehrer G, Meissner A, Bretschneider H J, Impedance Spectroscopy of tissue alterations during organ ischemia, 8th Inter. Conf. Elect. Bioimped., Kuopio, Finland, 176 - 178, 1992
7. Morucci J P, Aligne C, Chauveau N, Dumont P, Rigaud B and Cros S, Bioimpedance spectrometric variations of B16 tumour in mice during regression induced by dc current, IX Int. Conf. Elect. Bioimped., 219 - 222, 1995
8. Kruger W, Gersing E and Vaupel P, Electrical impedance spectroscopy for in vivo detection of structural changes in experimental tumours during local regression. IX Int. Conf. Elect. Bioimped., 225 - 228, 1995
9. Chauveau N., Dumont P, Aligne C, Rigaud B, Cross S and Morucci J P, In-vivo impedance spectrometry of MCF-7 tumours in nude mice: Measurements problems and examples. IX Int. Conf. Elect. Bioimped., 215 - 218, 1995
10. McAdams E T, and Jossinet J, Tissue impedance: a historical overview, *Physiol. Meas.*, **16**, A1 - A13, 1995.
11. Osypa M, and Gersing E, Tissue impedance spectra and the appropriate frequencies for EIT, *Physiol. Meas.*, **16**, A49 - A55, 1995
12. Ollmar S, Factors influencing the electrical properties of skin. *Innov. Tech. Biol. Med.*, **16**, 2, 137-142, 1995
13. Nicander I, Nyren M, Emtestam L, Ollmar S, Baseline electrical impedance measurements at various skin sites related to age and sex. *Skin Res & Tech*, **3**, 252-258, 1997
14. Ollmar S, Emtestam L, Electrical impedance applied to non-invasive detection of irritation in skin, *Contact Dermatitis*, **27**, 37-42, 1992
15. Ollmar S, Eek A, Sundstrom F, Emtestam L, Electrical impedance for estimation of irritation in oral mucosa and skin. *Medical Progress through Technology*, **21**, 29-37, 1995
16. Ollmar S, Nyren M, Nicander I and Emtestam L., Electrical impedance compared with other non-invasive bioengineering techniques and visual scoring for detection of irritation in human skin, *Brit Jour. Dermatology.*, **130**, 29-36, 1994
17. Schmukler R., Measurements of electrical Impedance of Living Cells in the Frequency domain., *Charge and Field Effects in Biosystems -2*, eds M. J. Allen, S F Cheary and M. F. Hawkrige, Plenum Press, 1989.
18. Freywald K H, Pliquett F, Schoberlein L, Pliquett U, Passive electrical properties of meat as a characterisation of it's quality, IX Int. Conf. Elect. Bioimped., 366 - 369, 1995.
19. Schmukler R, Impedance spectroscopy: The measurement of electrical impedance of biologic materials; *Electrical Trauma.*, eds. Lee R E, Cravalho E G, Burke J F, Cambridge Univ. Press, 239 - 253, 1992.
20. Schmukler R., Electrical Impedance of living cells: a modified four electrode approach. *IEEE, Engineering in Med. and Biol. Soc. New Orleans*, 1988

**For further details**

**Web:** <http://www.solartron.com/lap>

**E-mail:** [lab\\_info@solartron.com](mailto:lab_info@solartron.com)

**UK**

Victoria Rd., Farnborough  
Hampshire GU14 7PW England  
Telephone +44 (0) 1252 376666  
Fax +44 (0) 1252 544981

**USA**

964 Marcon Blvd. Suite 200  
Allentown, PA 18103, USA  
Telephone: +1 610-264-5034  
Fax +1 610-246-5329  
Toll-free 1-800 CALL SOL

**France**

37 rue du Saule Trapu  
91882 MASSY, Cedex, France  
Telephone +33 (0) 1 69 53 63 53  
Fax +33 (0) 1 60 13 37 06

**China**

Beijing Liaison Office  
Room 327, Ya Mao Building  
No. 16 Bei Tu Chen Xi Road  
Beijing 100101  
Peoples Republic of China  
Tel: +86 10-62381199 ext 2327  
Fax: +86 10-62384687